

EP-TRACER/38



EP System with USB connection.




The EP-TRACER/38 is a flexible PC controlled EP system, intended for all EP and ablation procedures in the clinical EP lab. In combination with a Laptop it can be used as a portable EP lab. It is the new generation of the Maastricht EP system. The EP-TRACER offers state of the art technology for the acquisition of 38 channels.

HARDWARE FEATURES


- 12 ECG channels
- **20 Intracardiac channels**
- 6 Auxiliary (pressure) channels
- Integrated 2 channel stimulator
- Stimulation on any intracardiac channel
- Optically isolated patient connections

SOFTWARE FEATURES

- Windows 2000/XP
- Clear and easy user interface
- Split screen with triggering on last stimulus
- Beat to beat triggered display
- Pace mapping
- Continuous storage on harddisk
- Event log
- Fast LaserJet printing

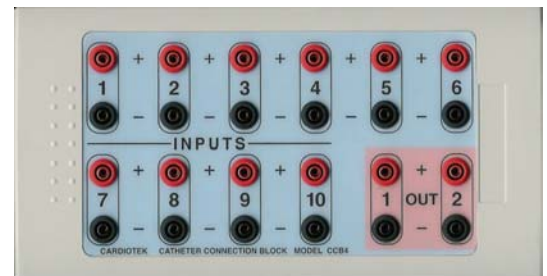

CCB4 - Catheter Connection Block

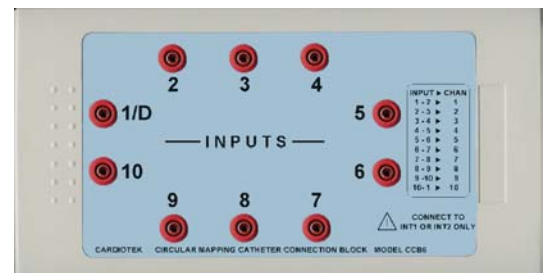

CCB6 - For circular mapping catheters

TECHNICAL SUMMARY

- Dynamic range: 20 bits
- Sample rate: 1 kHz / channel
- USB 1.1 and 2.0 compliant
- AC adapter 90-240VAC 50/60Hz
- Dimensions: 280 x 60 x 270 mm (W x H x D)
- Weight: 2,2 kg

CardioTek
Amerikalaan 70
6199AE Maastricht – Airport
The Netherlands

Tel: +31-43-3656006
Fax: +31-43-3656007
info@cardiotek.com
www.cardiotek.com